

Oxy-med[®]
O₂ for healthy living



Oxymed P2

Portable Oxygen Concentrator




Pulse Flow

1 to 5


Weight

1.98kg


Battery life

Up to
10 hrs
(double battery)

Effortless Portability

Oxy-med[®]
O₂ for healthy living


Enjoy the Freedom of Portable Oxygen Therapy, enabling you to travel, shop, and engage in activities with ease and peace of mind.

Weight
1.98 kg
only



Easy Control system with buttons and intuitive design

P2 offers effortless navigation with user-friendly buttons, putting ultimate relaxation control right at your fingertips

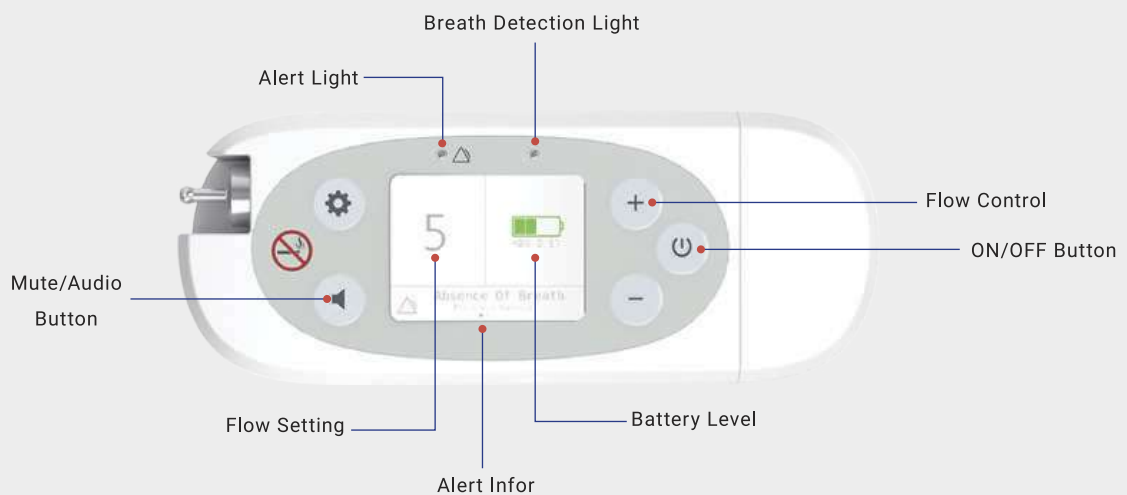

Approved for In-flight use

Drive with Oxygen

Keep your life-saving oxygen concentrator charged while driving, ensuring a constant supply of breathable air



Impressive Battery Backup up to 10 hours

Never worry about running out of power with our exceptional battery backup, giving you peace of mind for up to 10 hours of usage.

Pulse Setting	Duration (Single Battery)	Duration (Double Battery)
1	5 h 00 min	10 h 00 min
2	3 h 50 min	7 h 40 min
3	3 h 00 min	6 h 00 min
4	2 h 00 min	4 h 00 min
5	1 h 40 min	3 h 20 min

Recharge Time: 4 hours per battery

Uninterrupted Power

Use one, Charge the other

External Fast Charger*

Conveniently charge a backup battery while the device is in use.

Charging Time: 2.5 hours

*Additional purchase.



Specifications

Weight (with single battery)	1.98 kg
Size	8.70"L x 6.30"W x 3.35"H (includes battery)
Noise Level	49 dB @ Setting 2
Use	Designed for Continuous Use
Flow settings & Pulse Volume	
Oxygen Flow Settings (Pulse)	Pulse Dose Delivery System Upto 5
Maximum Minute Volume Output	Setting 1 - 210 mL / min Setting 2 - 420 mL / min Setting 3 - 630 mL / min Setting 4 - 840 mL / min Setting 5 - 1000 mL / min
Oxygen Concentration	90% - 3%/+6% at all settings
Warm up Time	2 minutes
Display Screen	2.8 inch
Battery & Power	
Battery Backup	Single Battery (Included) Upto 5 hours
	Double Battery (Additional Purchase) Upto 10 hours
Recharge Time (Per Battery)	4 hours 2.5 hours with external charger*
Battery Specs	Battery: Voltage: 14.4VDC Rated capacity: 6.7Ah

Power Supply	AC Input: 100 to 240VAC 50 to 60 Hz DC Input: 11-16 VDC DC Output: 19V 6.3A
Warranty & Quality Standards	
Quality Standards	FDA, CE, FAA Approved
In-Flight Use	Approved
Warranty	3 years (Main unit) 1 year (Battery and sieve beds)
What's in the Box	Oxymed P2 unit AC Adapter Car Charger Nasal Cannula Single Battery Sling bag Screwdriver 1 set of filters User manual Double battery is an additional purchase

*The device comes with single battery and additional purchase is required for double battery and external fast charger

Registered Office:

A brand by Medequip Healthcare Solutions LLP
Plot No: A-12, A-13, Apparel Park
Phase-1, Doddaballapur Industrial Area
Bengaluru - 561203, Karnataka

Branch Offices:

Delhi, Telangana, Odisha, Gujarat, Maharashtra,
Uttar Pradesh, Rajasthan, Madhya Pradesh Punjab,
Tamil Nadu, West Bengal, Kerala, Jharkhand,
Chhattisgarh and Bihar